

Qdot® Anti-Human CD3 Antibody Conjugates

Product	Form	Volume	Qdot®	Tests	Peak		Recommended
			Nanocrystal Conc.		Excitation (nm)	Emission (nm)	
Q10012	Qdot® 655	0.1 mL	1 µM	100 min.	405 (488)*	655	655/20
Q10484	Qdot® 655	0.025 mL	1 µM	25 min.	405 (488)*	655	655/20
Isotype Control: Mouse IgG2a							
Q10015	Qdot® 655	0.1 mL	1 µM	100 min.	405 (488)*	655	655/20

*Qdot® nanocrystals excite optimally in the UV to 405 nm range, but can also be excited with wavelengths shorter than their emission maximum, such as with a 488 nm laser

Product Description

Mouse monoclonal antibody to the human CD3 antigen

Clone: S4.1 (also known as 7D6)

Isotype: Mouse IgG2a

Lot No.: See label **Expiration:** See label

Buffer: 50 mM borate, 1 M betaine, pH 8.3

Preservative: 0.05% sodium azide. Sodium azide is an extremely toxic and dangerous compound particularly when combined with acids or metals. Properly dispose of solutions containing sodium azide.

Storage and Handling

Store reagents at 2–8°C. **Do not freeze.** Because Qdot® nanocrystals are conjugated to biological materials, some loss of activity may be observed with prolonged storage.

Qdot® nanocrystals are photostable, and do not need to be protected from light. However, if using Qdot® conjugates in combination with conventional fluorochrome conjugated antibodies, minimize light exposure during handling, incubation with cells, and prior to analysis. We recommend analysis of cells within 18 hours of staining. If dilute reagent is used, dilute only the quantity of reagent to be used within one day.

The Qdot® conjugate contains cadmium and selenium in an inorganic crystalline form. Dispose of the material in compliance with all applicable local, state, and federal regulations for disposal of these classes of material. For more information on the composition of these materials, consult the Material Safety Data Sheet.

Qdot® Nanocrystals

Qdot® nanocrystals are nanometer-scale atom clusters of semiconductor material which exhibit narrow and symmetrical emission bandwidths with very long Stokes shifts. The nanocrystals can be excited with any wavelength below their emission maximum, but are best excited by UV or violet light. Qdot® nanocrystals demonstrate an intrinsic brightness and photostability that can be many times greater than observed with other classes of fluorophores. These advantages make Qdot® nanocrystals powerful tools for antibody staining.^{1,2}

Product Characterization

Antigen Specificity: According to the literature this antibody recognizes the CD3 antigen.^{3–12} CD3 is comprised of 5 individual proteins that are responsible for T cell receptor expression and signaling. CD3 is expressed on thymocytes, T cells and some NK cells.

Leukocyte Workshop Status: Leukocyte Typing V

Product Use

Staining: Stain cells in any standard staining buffer, such as phosphate-buffered saline (PBS) with 1% bovine serum albumin (BSA). Use 1 µL of Qdot® antibody conjugate per 1×10^6 cells in a 100 µL staining volume (10 nM final concentration of Qdot® nanocrystals). Qdot® nanocrystal conjugates may be mixed with other antibodies, but use the diluted conjugates on the day of dilution. Because conditions may vary, you may need to determine the optimal amount of antibody used for each application. Qdot® nanocrystal conjugates can be used for surface staining applications with most conventional sample preparation reagents, such as Caltag™ Cal-Lyse™ and Caltag™ FIX & PERM® reagents, with minimal affect on fluorescence. We have observed that some batches of BD FACS™ Lysing Solution have been found to interfere with Qdot® nanocrystal fluorescence.


Instrument setup: Qdot® nanocrystals are excited optimally with UV or 405 nm light, although excitation can be obtained with any wavelength below the emission wavelength of a given nanocrystal. Make sure the cytometer has an appropriate emission filter for the Qdot® nanocrystal being used. The table above has filter recommendations; alternate filters can be used as long as they capture the emission maximum, but filter width impacts spectral overlap corrections.

Note: Qdot® nanocrystals can be used on cytometers that do not have UV or violet excitation sources as long as they have appropriate emission filters. Be sure to check for Qdot® nanocrystal emission in any channel that can capture nanocrystal emission off of other lasers on the cytometer.

Product Quality Control

Each lot has been tested by flow cytometry using human peripheral blood leukocytes (PBLs). This testing was performed using 1 µL of antibody per 1×10^6 cells in a 100 µL staining volume (10 nM final concentration of Qdot® nanocrystals). Qdot® nanocrystal concentration is assigned based on optical density. See reverse for representative flow cytometry data.

Qdot® Anti-Human CD3 Antibody Conjugates

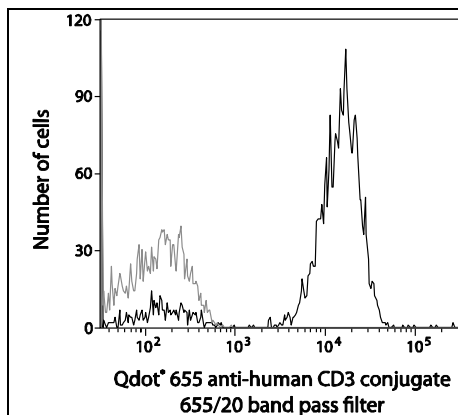

Log fluorescence intensity profile of human PBLs gated on lymphocytes and analyzed on a BD LSR II flow cytometer (BD Biosciences, San Jose, CA) using 405 nm excitation and the specified emission filter. The data file was analyzed using FlowJo™ software (Treestar, Inc., www.flowjo.com). The black line represents cells stained with anti-human CD3 antibody conjugate and the gray line represents unstained cells.

Note: Flow cytometric data shown may not necessarily have been generated using the enclosed lot of reagent. For this reason, and due to differences in flow cytometers and cytometer settings, results may vary from those illustrated above. We recommend titrating reagents to determine optimal conditions for use in your systems.

References

1. Telford, W. G. 2004. *Cytometry Part A* 61A:9.
2. Perfetto, S. P., P. K. Chattopadhyay, and M. Roederer. 2004. *Nature Reviews - Immunology* 4: 648.
3. Schlossman, S. F., L. Boumsell, W. Gilks, J. M. Harlan, T. Kishimoto, C. Morimoto, J. Ritz, S. Shaw, R. Silverstein, T. Springer, T. F. Tedder, and R. F. Todd eds. 1995. *Leukocyte Typing V*. Oxford University Press Inc., New York.
4. Reinherz, E. L., S. C. Meuer, K. A. Fitzgerald, R. E. Hussey, H. Levine, and S. F. Schlossman. 1982. *Cell* 30: 715.
5. Borst, J., S. Alexander, J. Elder, and C. Terhorst. 1983. *J. Biol. Chem.* 258: 5135.
6. Meuer, S. C., O. Acuto, K. A. Fitzgerald, R. E. Hussey, J. C. Hodgdon, S. F. Schlossman, and E. L. Reinherz. 1983. *Nature* 303: 808.
7. Meuer, S. C., K. A. Fitzgerald, R. E. Hussey, J. C. Hodgdon, S. F. Schlossman, and E. L. Reinherz. 1983. *J. Exp. Med.* 157: 705.
8. Van den Elsen, P., B. Shepley, J. Borst, J. E. Coligan, A. F. Markham, S. Orkin, and C. Terhorst. 1984. *Nature* 312: 413.
9. Furley, A. J., S. Mizutani, K. Weilbaecher, H. S. Dhaliwal, A. M. Ford, L. C. Chan, H. V. Molgaard, B. Toyonaga, T. Mak, P. van den Elsen, D. Gold, C. Terhorst, and M. F. Greaves. 1986. *Cell* 46: 75.
10. Gold, D. P., J. M. Puck, C. L. Petty, M. Cho, J. Coligan, J. Woody, and C. Terhorst. 1986. *Nature* 321: 4.
11. Oettgen, H. C., and C. Terhorst. 1987. *Hum. Immunol.* 18: 187.
12. Clevers, H., B. Alarcon, T. Wileman, and C. Terhorst. 1988. *Ann. Rev. Immunol.* 6: 629.

Related Products

Catalog no.	Product Name	Unit Size
GAS-010	Cal-Lyse™ Whole Blood Lysing Solution	25 mL
GAS-010S-100	Cal-Lyse™ Whole Blood Lysing Solution	100 mL
HYL-250	High-Yield Lyse Fixative	500 mL
GAS001S-5	FIX & PERM® Reagent A (Individual)	5 mL
GAS001S-100	FIX & PERM® Reagent A (Bulk)	100 mL
GAS002S-5	FIX & PERM® Reagent B (Individual)	5 mL
GAS002S-100	FIX & PERM® Reagent B (Bulk)	100 mL
GAS-003	FIX & PERM® Reagents	50 tests

Limited Use Label License: Research Use Only

The purchase of this product conveys to the purchaser the limited, non-transferable right to use the purchased amount of the product only to perform internal research for the sole benefit of the purchaser. No right to resell this product or any of its components is conveyed expressly, by implication, or by estoppel. This product is for internal research purposes only and is not for use in commercial services of any kind, including, without limitation, reporting the results of purchaser's activities for a fee or other form of consideration. For information on obtaining additional rights, please contact outlicensing@lifetech.com or Out Licensing, Life Technologies, 5791 Van Allen Way, Carlsbad, California 92008.

©2011 Life Technologies Corporation. All rights reserved. This information is subject to change without notice.

The trademarks mentioned herein are the property of Life Technologies Corporation or their respective owners.

FACS is a trademark of BD Biosciences. FlowJo is a trademark of Treestar, Inc.

For support visit www.lifetechnologies.com/support or email techsupport@lifetech.com
www.lifetechnologies.com